

WHITE BLOOD CELL MIGRATION STUDY (In-111-WBCs, Tc-99m-HMPAO-WBCs)

Overview

- The White Blood Cell Migration Study demonstrates the distribution of labeled autologous white blood cells within the body at various times after intravenous injection. When labeled with technetium-99m, the study can be completed in several hours and is optimal for imaging inflammatory bowel disease. When labeled with indium-111, the study can be extended over two days and is preferred for the diagnosis of vertebral body osteomyelitis. White blood cells labeled with either radioisotope may be used for imaging osteomyelitis outside of the spine and soft tissue abscesses.

Indications

- Detection of abscesses and infection in soft tissues, particularly in patients without localizing findings (1-3).
- Diagnosis of infection in the skeleton (4-8).
- Evaluation of inflammatory bowel disease (9-11).

Examination Time

- Initially: 15 minutes for withdrawal of 50 mL of blood.
- At 1 hour: 15 minutes for injection of the radiopharmaceutical.
- Later:
 - > In-111-WBCs: 1 hour for imaging at 6 and/or 24 hours (12). (Delayed images beyond the routine set of images may be needed.)
 - > Tc-99m-HMPAO-WBCs: 1 hour for imaging at 2 hours (3,13).

Patient Preparation

- None (antibiotics do not appear to interfere) (3).

Equipment & Energy Windows

- Gamma camera: Large field of view, preferably with dual heads.
- Collimator:
 - > In-111-WBCs: Medium energy, parallel hole.
 - > Tc-99m-HMPAO-WBCs: Low energy, high resolution, parallel hole.

- Energy windows:
 - > In-111-WBCs:
 - θ One pulse height analyzer: 156 to 272 keV.
 - θ Two pulse height analyzers: 20% windows centered at 171 and 245 keV.
 - > Tc-99m-HMPAO-WBCs: 20% window centered at 140 keV.

Radiopharmaceutical, Dose, & Technique of Administration

- Radiopharmaceutical (1,3,13-15):
 - > In-111-white blood cells - radiopharmaceutical of choice for infection of the spine, equal for infection of rest of skeleton and soft tissue infection.
 - > Tc-99m-HMPAO-WBCs - radiopharmaceutical of choice for inflammatory bowel disease, equal for soft tissue infection (1,3,13).
- Dose:
 - > In-111-WBCs: 500 μ Ci (18.5 MBq) (8).
 - > Tc-99m-HMPAO-WBCs: 10-25 mCi (370 MBq) (3).
- Technique of administration: Standard intravenous injection.

Patient Position & Imaging Field


- Patient position: Supine.
- Imaging field: Usually entire torso; check with the nuclear medicine

Acquisition Protocol

- Initial imaging:
 - > In-111-WBCs : 24 hours (occasionally 6 hours).
 - > Tc-99m-HMPAO-WBCs: 2 hours for (3).
- ANT and POST images are acquired of the torso, and of the extremities as indicated; occasionally other projections are obtained.
 - θ For moving acquisition: Use a camera/table motion of approximately 5-10 cm/min.
 - θ Static acquisition protocol: Acquire images for approximately 5 minutes each.
- SPECT imaging - use routinely for limited studies and areas of uncertainty (16):
 1. Image acquisition parameters:
 - a) degrees of rotation: 360°.
 - b) number of images: 64.
 - c) time per image: 20 seconds.

- Delayed imaging (may be needed to increase the certainty of diagnosis (1)):
 - > In-111-WBCs: 48 hours.
 - > Tc-99m-HMPAO-WBCs: 24 hours for Tc-99m-HMPAO-WBCs (16),

Protocol Summary Diagrams

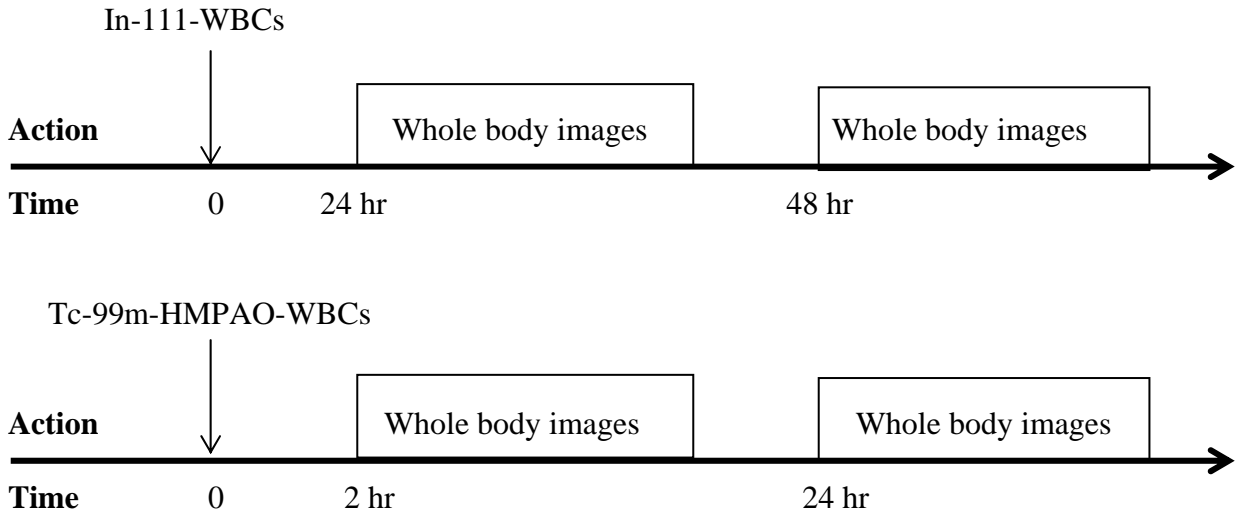

Data Processing

- SPECT image reconstruction:
 1. The exact procedure for processing SPECT images depends on the computer software being used. This varies with the manufacturer and, in general, the manufacturer's protocol should be followed.
 2. The reconstruction process in general terms is:
 - a) correct the 64 planar images for uniformity (camera non-uniformity) using a high count, e.g. 30 million count, cobalt-57 flood acquisition.
 - b) check the images for patient motion and apply a motion correction algorithm if indicated and if available.
 - c) if the entire field of view is not of interest, indicate the region that is of interest so that computer time is not expended reconstructing tomograms outside the region of interest.
 - d) specify the filters to be used in the reconstruction process and the pixel thickness of the tomogram (usually 1 or 2 pixels).
 - e) reconstruct transverse, sagittal, and coronal image.

Optional Maneuvers

- Bone marrow imaging with Tc-99m-sulfur colloid: May be performed in conjunction with In-111-WBC studies for osteomyelitis to increase

- Pethidine to reduce hepatobiliary excretion: When imaging inflammatory bowel disease with Tc-99m-HMPAO-WBCs, can give pethidine (an opiate) to constrict the sphincter of Oddi and reduce excretion of Tc-99m-WBC by products (16).
- Bone mineral imaging with Tc-99m-HMP or MDP: May be performed in conjunction with In-111-WBC studies for osteomyelitis to increase

Principle Radiation Emission Data - In-111 (21)

- Physical half-life = 2.83 days.

Radiation	Mean % per disintegration	Mean energy (keV)
Gamma-2	90.2	171.3
Gamma-3	94.0	245.3

Dosimetry - In-111-White Blood Cells (22,23)

Organ	rads/500 μ Ci	mGy/18.5 MBq
Spleen	11.9	119.0
Liver	1.9	19.0
Lungs	0.8	7.5
Marrow	0.3	3.3
Total body	0.2	2.0
Ovaries	0.1	0.7
Testes	0.01	0.1

Principle Radiation Emission Data - Tc-99m (24)

- Physical half-life = 6.01 hours.

Radiation	Mean % per disintegration	Mean energy (keV)
Gamma-2	89.07	140.5

Dosimetry - Tc-99m-HMPAO-White Blood Cells (25)

Organ	rads/10 mCi	mGy/370 MBq
Spleen	2.3	23.3
Liver	1.6	16.3
Kidneys	1.4	13.7
Bladder wall	0.6	6.3
Lungs	0.4	3.7
Marrow	0.1	1.5
Total body	0.1	1.5
Ovaries	0.1	0.7
Testes	0.04	0.4

References

1. Palestro CJ, Love C, Tronco GG, et al: Role of radionuclide imaging in the diagnosis of postoperative Infection. RadioGraphics 20:1649-1660, 2000.
2. Peters AM: The utility of Tc-99m HMPAO-leukocytes for imaging infection. 24:110-127, 1994.
3. Datz FL: Abdominal abscess detection: Gallium, In-111, and Tc-99m-labeled leukocytes, and polyclonal and monoclonal antibodies. Sem Nucl Med 26:51-64, 1996.
4. Krznicaric E, De Roo M, Verbruggen A, et al: Chronic osteomyelitis: Diagnosis with technetium-99m-d,l-hexamethylpropylene amine oxime labelled leucocytes. Eur J Nucl Med 23:792-797, 1996.
5. Devillers A, Moisan A, Jean S, et al: Technetium-99m hexamethylpropylene amine oxime leucocyte scintigraphy for the diagnosis of bone and joint infections: A retrospective study in 116 patients. Eur J Nucl Med 22:302-307, 1995.
6. Palestro CJ, Roumanas P, Sawyer AJ, et al: Diagnosis of musculoskeletal infection using combined In-111 labeled leukocyte and Tc-99m SC marrow imaging. Clin Nucl Med 17:269-273, 1992.
7. Palestro CJ, Sawyer AJ, Kim CK, et al: Infected knee prosthesis: Diagnosis with In-111 leukocyte, Tc-99m sulfur colloid, and Tc-99m MDP imaging. Radiology 179:645-648, 1991.
8. Love C, Marwin SE, Tomas MB, et al: Diagnosing infection in the failed joint replacement: A comparison of coincidence detection F-18-FDG and In-111-labeled leukocyte/Tc-99m-sulfur colloid marrow imaging. J Nucl Med 45:1864-1871, 2004.
9. Allan RA, Sladen GE, Bassingham S, et al: Comparison of simultaneous Tc-99m-HMPAO and In-111 oxine labelled white cell scans in the assessment of inflammatory bowel disease. Eur J Nucl Med 20:195-200, 1993.
10. Bennink RJ, Peeters M, Rutgeerts, et al: Evaluation of early treatment response and predicting the need for colectomy in active ulcerative colitis with Tc-99m-HMPAO white blood cell scintigraphy. J Nucl Med 45:1698-1704, 2004.
11. Charron M, del Rosario FJ, Kocoshis SA: Pediatric inflammatory bowel disease: Assessment with scintigraphy with Tc-99m white blood cells. Radiology 212:507-513, 1999.
12. Arndt JW, van der Sluis Veer A, Blok D, et al: Prospective comparative study of technetium-99m-WBCs and indium-111-granulocytes for the examination of patients with inflammatory bowel disease. J Nucl Med 34:1052-1057, 1993.
13. Mountford PJ, Kettle AG, O'Doherty MJ, et al: Comparison of technetium-99m-HM-PAO leukocytes with indium-111-oxine leukocytes for localizing intraabdominal sepsis. J Nucl Med 331:311-315, 1990.
14. Schauwecker DS, Burt RW, Park HM, et al: Comparison of purified indium-111 granulocytes and indium-111 mixed leukocytes for imaging of infections. J Nucl Med 29:23-25, 1988.

15. Bertrand-Caix J, Freyburger G, Bordenave L, et al: Functional upregulation of granulocytes labeled with technetium-99m-HMPAO and indium-111-oxinate. J Nucl Med 37:863-868, 1996.
16. Larikka MJ, Ahonen AK, Junila JA, et al: Extended combined tc-99m-white blood cell and bone imaging improves the diagnostic accuracy in the detection of hip replacement infections. Eur J Nucl Med 28:288-293, 2001.
17. Palestro CJ, Kim CK, Sawyer AJ, et al: Total-hip arthroplasty: Periprosthetic indium-111-labeled leukocyte activity and complementary technetium-99m-sulfur colloid imaging in suspected infection. J Nucl Med 31:1950-1955, 1990.
18. Seabold JE, Nepola JV, Marsh JL, et al: Postoperative bone marrow alterations: Potential pitfalls in the diagnosis of osteomyelitis with In-111-labeled leukocyte scintigraphy. Radiology 180:741-747, 1991.
19. Achong DM, Oates E: The computer-generated bone marrow subtraction image: A valuable adjunct to combined In-111 WBC/Tc-99m in sulfur colloid scintigraphy for musculoskeletal infection. Clin Nucl Med 19:188-193, 1994.
20. Johnson JA, Christie MJ, Sandler MP, et al: Detection of occult infection following total joint arthroplasty using sequential technetium-99m HDP bone scintigraphy and indium-111 WBC imaging. J Nucl Med 29:1347-1353, 1988.
21. 49-In-111: In MIRD: Radionuclide Data and Decay Schemes, DA Weber, KF Eckerman, AT Dillman, JC Ryman, eds, Society of Nuclear Medicine, New York, 1989, pp 196-197.
22. Williams LE, Forstrom LA, Weiblen BJ, et al: Human dosimetry of In-111 granulocytes. In In-111 Labeled Platelets and Leukocytes. Proceedings of the 2nd Annual Symposium, HW Wahner, DA Goodwin, eds, Rochester, MN, Mayo Clinic, October, 1981, pp 175-188.
23. Gainey MA, Siegel JA, Smergel AM, et al: Indium-111-labeled white blood cells: Dosimetry in children. J Nucl Med 29:689-694, 1988.
24. 43-Tc-99m: In MIRD: Radionuclide Data and Decay Schemes, DA Weber, KF Eckerman, AT Dillman, JC Ryman, eds, Society of Nuclear Medicine, New York, 1989, pp 178-179.
25. Skretting A, Benestad HB, Sundrehagen E: Whole-body distribution of Tc-99m-labelled autologous human granulocytes and radiation dose to cells and organs. Eur J Nucl Med 14:1-7, 1988.

Normal Findings

- > Guze BH, Webber MM, Hawkins RA, et al: Indium-111 white blood cell scans: Sensitivity, specificity, accuracy, and normal patterns of distribution. Clin Nucl Med 15:8-10, 1990.
- > Oswald SG, Van Nostrand D, Sarvory CG, et al: The acetabulum: A prospective study of three-phase bone and indium white blood cell scintigraphy following porous-coated hip arthroplasty. J Nucl Med 15:274-280, 1990.
- > Datz FL: Abdominal abscess detection: Gallium, In-111, and Tc-99m-labeled leukocytes, and polyclonal and monoclonal antibodies. Sem Nucl Med 26:51-64, 1996.

- > Kipper SL: Radiolabeled leukocyte imaging of the abdomen: Nuclear Medicine Annual pp 81-128, 1995.